



Incavo Tibial Component Revision Osteotomes

Designed by Stephen J. Incavo, MD

Designed to help break the posterior cement-bone interface when removing a cemented tibial TKA component

Also used to help break the posterior implant-bone interface when removing a cementless tibial TKA component.



PRODUCT NO'S:

3621-01 [Standard Osteotome]


Blade Length: 10mm
Blade Width: 1/2" (12.7mm)
Blade Offset: 3/4" (19.1mm)
Overall Length: 8.5"

3621-02 [Medium Osteotome]

Blade Length: 14mm
Blade Width: 1/2" (12.7mm)
Blade Offset: 3/4" (19.1mm)
Overall Length: 8.5"

3040 [Slap Hammer]

1015 [Sterilization Case]

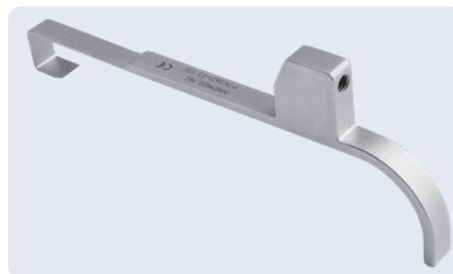

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

www.innomed.net



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

INNOMED